

SUPRENO®


POWDER-FREE NITRILE EXAMINATION GLOVES

• STANDARD EXAMINATION GLOVE



Features & Benefits:

- Soft and durable nitrile formulation for improved elasticity and greater comfort
- Interior polymer coating makes glove easy to don
- 100% nitrile to reduce the potential for adverse reactions associated with latex protein
- Textured fingers allow for reliable grip in wet and dry conditions
- Gradient dipping provides added protection where you need it most
- Meets all current FDA and EN/ISO standards assuring consistent quality in every box, every time
- EN 374-2, EN 374-3 certified

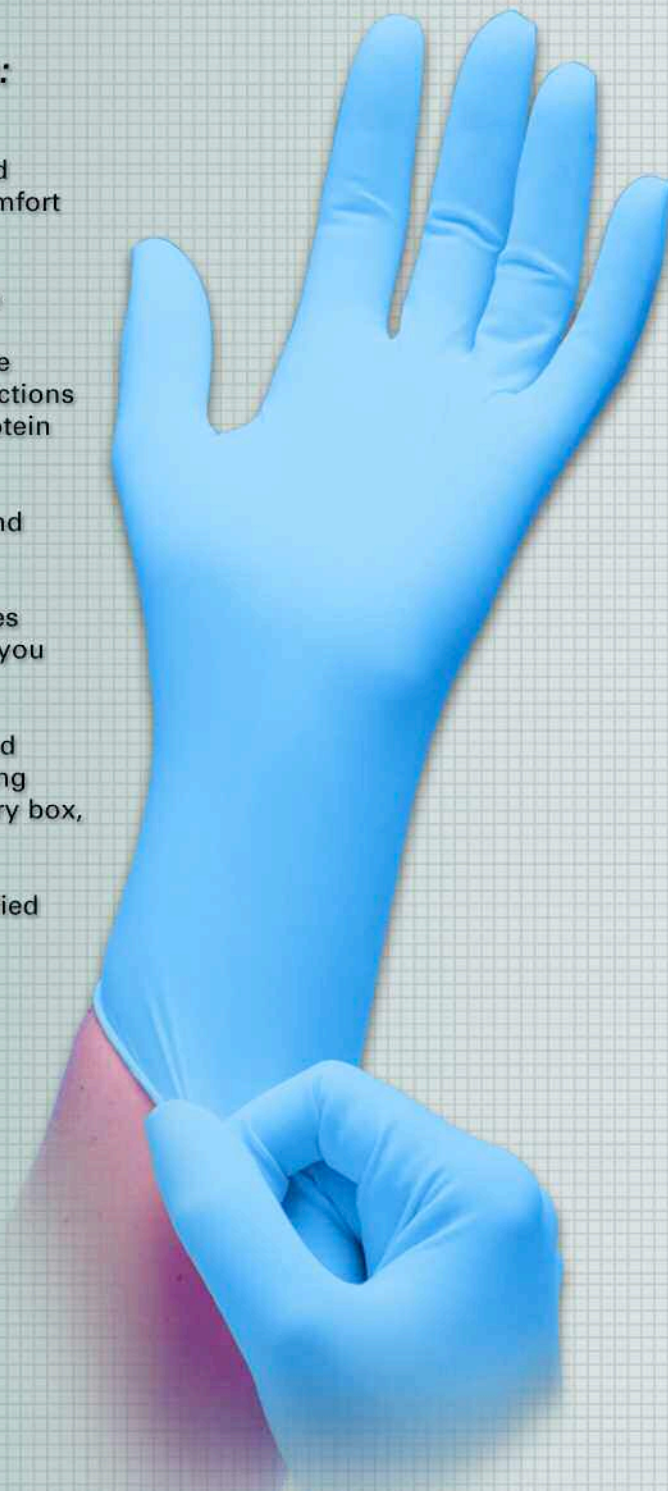

Exceptional comfort in a tough nitrile glove.

Supreno® is a soft, comfortable nitrile examination glove formulated for superior strength and durability. More comfortable than nitriles of the past, true-fitting **Supreno®** won't fatigue hands over long periods of wear. Enhanced tactile sensitivity and textured fingers offer better gripping in wet and dry conditions. Interior polymer coating makes **Supreno®** easy to don. Online chlorination processing means more consistent quality in every box, every time.

EN-374



MICROFLEX®
THE MOST TRUSTED NAME IN GLOVES®

www.microflex.com/int

REGISTERED
ISO 13485

SUPRENO®

POWDER-FREE NITRILE EXAMINATION GLOVES

STANDARD EXAMINATION GLOVE



GLOVE PROFILE

Supreno*

Length	245mm / 9.6 in.
Cuff Thickness	4.3 mils
Palm Thickness	5.5 mils
Finger Thickness	7.9 mils
Force at Break	min. 9
Elasticity/Elongation (%)	min. 500
Powder Mass	max. 2 mg / glove

GLOVE ATTRIBUTES

Material:	Nitrile; not made from natural rubber latex
Type:	Non-sterile
Shape:	Ambidextrous
Color:	Blue
Sizes:	Extra-Small, Small, Medium, Large, Extra-Large
Cuff:	Beaded
Glove Interior:	Polymer-coated, powder-free
Glove Exterior:	Textured fingers
Usage:	For single use only
Packaging:	100 gloves per dispenser box, by weight; 10 boxes per case

* ISO11193-1; EN455 I, II, III; EN 374-2, EN 374-3; CE Class I

Powder-free synthetic examination gloves

See for yourself why **Microflex®** gloves are the first choice of exam glove users nationwide. For samples or more information on any of the fine latex and synthetic products from **Microflex®**, contact your **Microflex®** distributor or for the name of a distributor in your area, call:

775-746-6600 ext. 1726


Or visit our website at:

www.microflex.com/int

Fax: 775-746-6586

PRODUCT CODE	RANGE OF SIZE	DESCRIPTION	CASE QUANTITY
SU-INT-XS	EXTRA-SMALL	Nitrile, Non-sterile, Examination	1000
SU-INT-S	SMALL	Nitrile, Non-sterile, Examination	1000
SU-INT-M	MEDIUM	Nitrile, Non-sterile, Examination	1000
SU-INT-L	LARGE	Nitrile, Non-sterile, Examination	1000
SU-INT-XL	EXTRA-LARGE	Nitrile, Non-sterile, Examination	1000

Caution Statement: Components used in making these gloves may cause allergic reactions in some users. Follow your institution's policies for use.

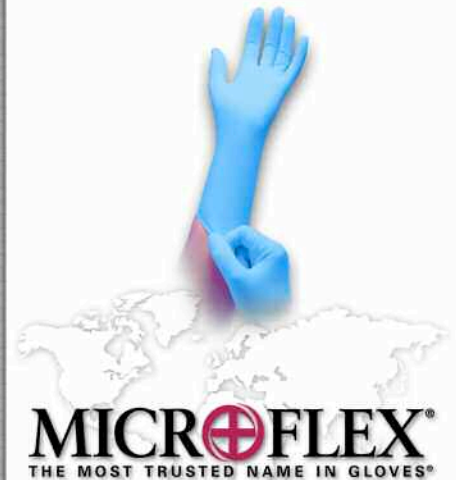

MICROFLEX®
THE MOST TRUSTED NAME IN GLOVES®

www.microflex.com/int

©2007 Microflex Corporation
P.O.Box 32000
Reno, NV 89533-2000
All Rights Reserved.
07.00003 AL 061907

EUROPEAN PATENT NO. 0456333.